



TissueScope™ LE120 Slide Scanner

Best-in-class slide scanning for medium volume labs.

The TissueScope LE120 is designed for labs that require the best-in-class image quality of our versatile TissueScope LE scanner, but with higher capacity. Each of the scanner's 10 slide holders accommodates up to twelve 1" x 3" slides, allowing for continuous scanning of up to 120 brightfield slides at once.

Key Features

- Sharp, clear image quality up to 40X magnification
- Scan up to 120 brightfield slides continuously with automatic setup and non-stop operation
- Scan standard 1" x 3" and 2" x 3" slides as well as whole mount slides in any size up to 6" x 8"
- Includes HuronViewer software for image viewing
- Optional TissueView image server and management software platform
- Pair with the TissueSnap preview station for enhanced workflow optimization

Technical Specifications


Imaging Mode	Brightfield
Slide Formats	25mm x 75mm (1" x 3") 50mm x 75mm (2" x 3") Whole mount slides up to 150mm x 200mm (6" x 8")
Optical Resolution	0.75 NA
Resolution (µm/pixel)	0.25 at 40X (2.0 at 5X) 0.50 at 20X (5.0 at 2.5X) 1.00 at 10X (10.0 at 1.0X)
Total Time to View¹	<5 minutes per slide ²
File Format	Non-proprietary 24-bit RGB Pyramidal BigTIFF Uncompressed, JPEG or JPEG2000 Compression Export to JPEG, Flat TIFF and LZW Compressed TIFF
Dimensions & Weight	89.0cm (35.0") Width 60.5cm (23.8") Length 74.5cm (29.3") Height 180Kg (396lbs) Weight
Operating Temperature	10°C - 40°C (50°F - 104°F)
Certification	CSA
Warranty	One year

¹ Total Time to View includes slide loading, previewing, scanning, unloading, compression and file saving.
² 15mm x 15mm @ 20X magnification.

Mix and match slide formats – standard to whole mount

Our unique slide holder has the highest slide capacity and can accommodate slide formats from 1" x 3" to whole mount up to 6" x 8".

And with our secure clamping system, slides are securely seated in the slide holder. This keeps your precious specimens and fragile glass safe and gives you peace of mind.

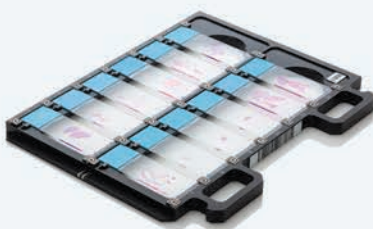

The optional TissueView digital image management system integrates with your LIS to make your workflow easier.



Designed with and by pathologists, TissueView provides an intuitive and customizable workflow. On a single integrated workstation, review cases and generate reports, saving time by accessing slides from anywhere. Annotate directly onto images and measurements as you view. Share cases easily for review and quality control and complete live collaborations with multiple remote attendees for instant consultation.

Optimize your workflow with the TissueSnap preview station

The TissueSnap preview station works with our TissueScope scanners to offload time consuming tasks such as preview scanning and batch pre-processing so the scanner can work at peak throughput and performance.



About Huron Digital Pathology

We believe that innovation is the key to making digital pathology a ubiquitous reality. But for us, innovation is about more than just delivering great image quality and fast scanning speeds. It's also about designing products that are simple to use, easy to integrate with your workflow, and have attractive price-performance characteristics.

Based in Waterloo, Ontario, Canada, Huron Digital Pathology has a 20 year history designing sophisticated imaging instrumentation. Our end-to-end digital whole slide scanning solutions for digital pathology incorporate our award-winning TissueScope™ digital slide scanners; TissueView™ image viewing, sharing and management platform; and our workflow-enhancing accessories, which include our innovative TissueSnap™ preview scanning station.

The TissueScope™ Scanner is intended for research and educational purposes only. TissueScope™, TissueView™, TissueSnap™ and Huron Digital Pathology™ are trademarks of Huron Digital Pathology. The TissueScope™ Scanner is based on one or more of the following U.S. Patents: 8,896,918 B2/8,655,043/7,218,446 B2 / 6,072,624 / 5,760,951 / 5,737,121 / 5,532,873 / 5,381,211. U.S. and International patents pending. Huron Digital Pathology is an operating division of Huron Technologies International Inc.

Distributed by



MEL: (03) 9480 4999
SYD: (02) 9705 8059
sales@scitech.com.au
www.scitech.com.au